



**elliquence**  
Less Is More®

# Trigger·Flex<sup>Mini</sup>®

## Exclusive Treatment for Endoscopic Spine Surgery

- Compatible with most Spine Scopes\*
- Safe and precise Soft Tissue/Nucleus Ablation and Coagulation
- Help provide a clear picture to determine anatomic landmarks
- Quick Control of bleeding to regain visualization
- Safe and precise Annulus Modulation
- Steerable/Navigational Access
- Tactile Feedback/Probing
- Radiopaque



\* Spine Scopes must have an inner working channel diameter of 2.4mm

## Product Specifications and Technology Overview



**Modes for Trigger-Flex<sup>®</sup> Mini with Surgi-Max<sup>®</sup>**  
**Bipolar Hemo**  
**Bipolar Turbo**



Used exclusively with Surgi-Max<sup>®</sup>  
 Radiowave Technology




Articulated Surgical Precision  
 for targeted areas



Low Heat Temperatures  
 allow minimal damage  
 to surrounding tissue



Mini Bipolar Tip  
 for focused energy

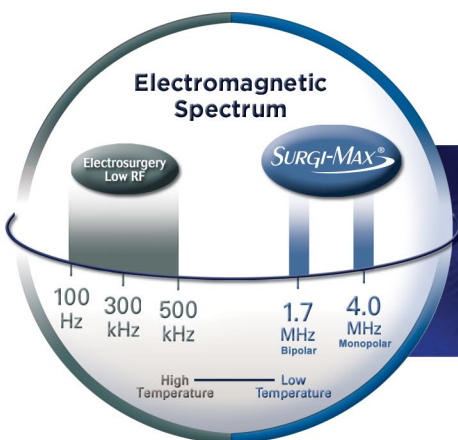


## Our Technology...See The elliquance<sup>®</sup> Difference

elliquance<sup>®</sup> has manufactured innovative radiofrequency systems and electrodes since 1969. The Trigger-Flex<sup>®</sup> handpiece combined with the Surgi-Max<sup>®</sup> High Frequency Radiowaves is designed for minimally invasive procedures such as endoscopic spine surgery. While laser and most electro-surgical devices are heat driven and imprecise, Surgi-Max<sup>®</sup> provides unparalleled surgical precision at a much lower temperature allowing it to gently target tissue with minimal damage to surrounding healthy tissue. A, B, C, D, E



**DTFM-31**

Working Length: 31cm  
 Outer Diameter: 2.3mm

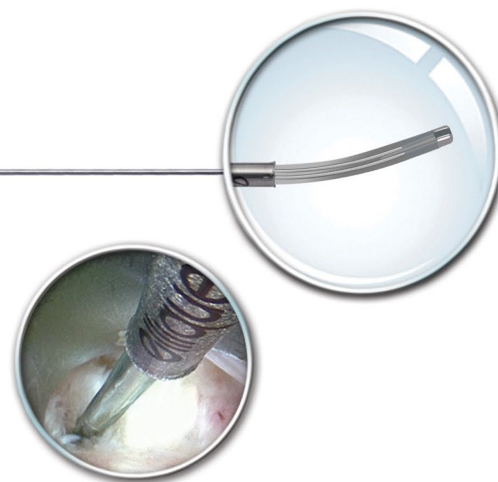

**Warning and precautions:**

- A) Hazardous electrical output. This equipment is for use only by qualified personnel.
- B) Read the instruction manual in detail before use.
- C) Use of accessories or cables other than elliquance specified, as replacement parts for internal components, may result in increased emissions or decreased immunity of the elliquance RF generator.
- D) The coupling of any non-elliquance manufactured product, with an elliquance manufactured product invalidates all elliquance product warranties and all elliquance legal liabilities associated with such use.
- E) Only elliquance supplied and/or approved accessories should be used to ensure proper operation. The elliquance generator and its accessories are identified with the manufacturer information and reference identification number for proper traceability.